


Separation of Clopidogrel Bisulfate and Related Substances on Ultron ES-OVM (L57)

(columns available from Hichrom Ltd)

Clopidogrel Bisulfate

Column	ULTRON ES- OVM 150 x 4.6 mm I.D.
Sample	Clopidogrel Bisulfate <i>RS</i> : 0.5µ g/ mL
	Related Compound A <i>RS</i> : 1µ g/ mL
	Related Compound B <i>RS</i> : 3µ g/ mL
	Related Compound C <i>RS</i> : 5µ g/ mL
Mobile Phase	10mM KH ₂ PO ₄ / Acetonitrile = 75 / 25
Flow Rate	1.0 mL / min
Injection Volume	10µ L
Temperature	25 degrees celsius
Detection	220nm

 SHINWA CHEMICAL INDUSTRIES, LTD.

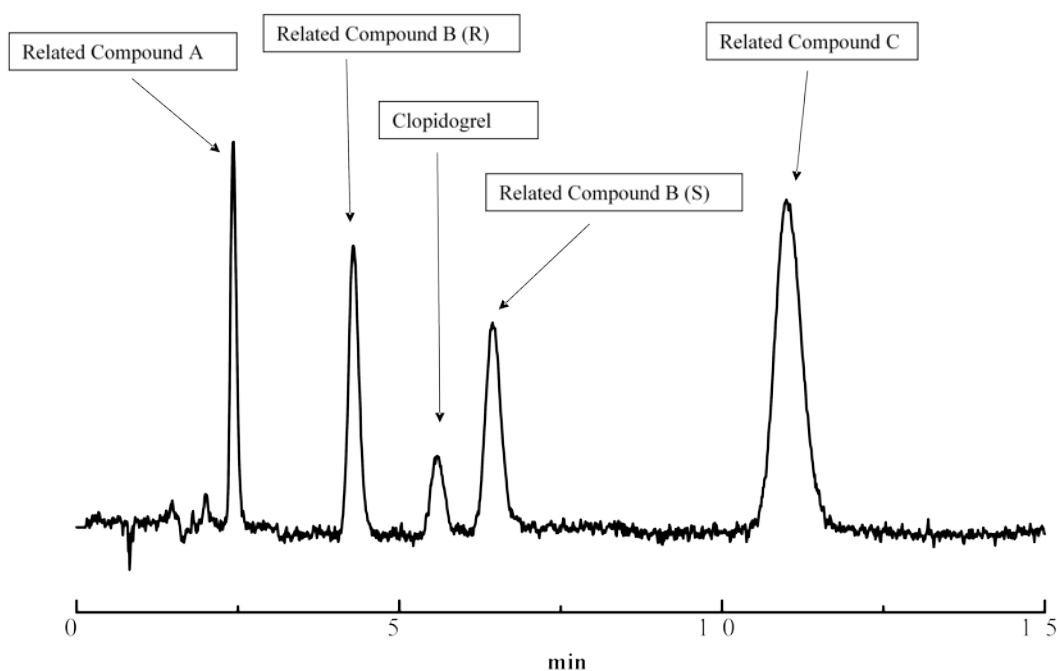

Column Details

Description	Ultron ES-OVM 5µm, 150 x 4.6mm
Part Number	ESOVM-1546

